


10mm Square HF Compatible Chip Holder Instruction



10mm Square Chip Holder is not compatible with Piranha Etch



Teflon Base Top View

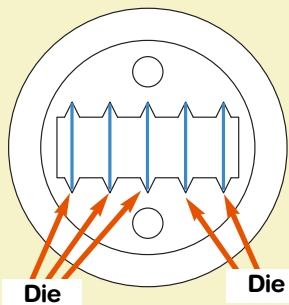

Teflon Base

Polypropylene Cover

#8768E 10mm Square HF Compatible Chip Holder with Chip Holder Handle

Twist to unscrew the Polypropylene Cover of #8768E

Top View

Die

Die

Load 10mm square chip into each slot of the Teflon Base. You may load up to 5 die into Teflon Base.

Front View

Twist to screw the Polypropylene Cover of #8768E

#8768E can be transported using handle provided or a tweezer